

The TCIA Plugin and Servlet API


Table of Contents

1	Installing the TCIAPlugin in a CTP instance:	3
2	The config.xml file	4
3	User Accounts	6
4	The TCIAServlet API	7
4.1	Updating the DicomAnonymizer Lookup Table	7
4.2	Submitting Files to the Import Pipeline	8
4.2.1	Listing Files on the Server	8
4.2.2	Submitting Files to the DirectoryImportService	8
4.2.3	Submitting Files to the DicomImportService	8
4.3	Viewing the Status of the Import Pipeline	8
4.4	Viewing the List of Files Ready for Anonymization	9
4.5	Submitting Files for Anonymization	9
4.6	Viewing the List of Files Ready for Export	10
4.7	Exporting Files	11
4.8	The Manifest	11
4.8.1	Clearing the Manifest	11
4.8.2	Getting the Manifest Status	11
4.8.3	Listing the Manifest with PHI as XML	11
4.8.4	Listing the Manifest without PHI as XML	12
4.8.5	Listing the Manifest with PHI as CSV	13
4.8.6	Listing the Manifest without PHI as CSV	13
4.8.7	Listing the Manifest with PHI as XLSX	13
4.8.8	Listing the Manifest without PHI as XLSX	13
4.8.9	Exporting the Manifest	13
4.9	Image Functions	13
4.9.1	Listing the Elements of a DICOM File	13
4.9.2	Getting a JPEG Image from a DICOM File	14
4.10	Special Functions	14

4.10.1	Getting the List of File System Roots on the Server	14
4.10.2	Getting the Available Space on the Server	14
4.10.3	Getting the Export Queue Size.....	14
4.10.4	Getting the URL of the Quarantine Manager Servlet	15
4.10.5	Getting a Summary of the DicomAnonymizer Quarantine.....	15
4.10.6	Shutting Down the Server	15
4.11	Testing Functions.....	15

1 Installing the TCIAPlugin in a CTP instance:

1. Install or upgrade to the latest CTP.
2. Put the TCIAPlugin.jar file in the **CTP/libraries** folder.
3. Put the special TCIA version of the config.xml file (see below) in the **CTP** folder.
4. Put any web pages, css files, Javascript files, and wizard files in the **CTP/ROOT** directory tree, being careful not to overwrite the index.xml file.
5. Run the Launcher program:

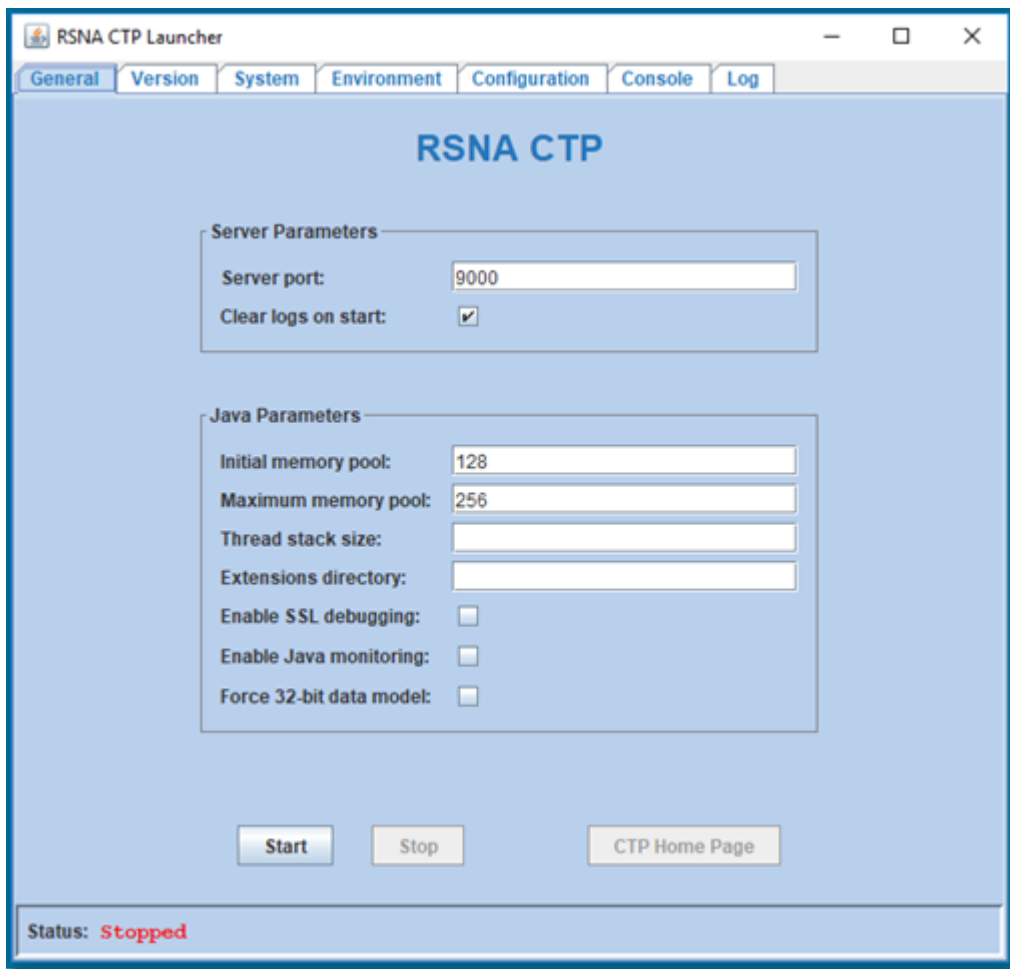

6. The config.xml file puts the CTP server on port 9000. It can be changed by altering the **Server port** field.
7. Click the **Start** button.
8. The **CTP Home Page** button will then be enabled. Click it.
9. The default browser will start and display the home page.
10. Log in as username=**admin**, password=**password**.
11. To shut down the server, click the **Stop** button on the Launcher window (or just close the Launcher).

2 The config.xml file

The config.xml for the TCIA project is:

```
<Configuration>
  <Server
    maxThreads="20"
    port="9000"/>
  <Plugin
    anonymizerID="CollectionDicomAnonymizer"
    anonymizerInputID="CollectionAnonymizerInput"
    anonymizerStorageID="CollectionAnonymizerStorage"
    class="edu.uams.tcia.TCIAPlugin"
    exportInputID="CollectionExportInput"
    exportManifestLogID="CollectionManifestLog"
    exportOutputID="CollectionExportOutput"
    id="Collection"
    importInputID="CollectionImportInput"
    importStorageID="CollectionImportStorage"
    name="CollectionTCIAPlugin"
    root="Collection/TCIAplugin"/>
  <Plugin
    class="edu.uams.tcia.ManifestLogPlugin"
    id="CollectionManifestLog"
    name="CollectionManifestLog"
    root="Collection/ManifestLog"/>
  <Pipeline
    name="CollectionImport"
    root="Collection/import">
    <DirectoryImportService
      class="org.rsna.ctp.stdstages.DirectoryImportService"
      id="CollectionImportInput"
      import="roots/DirectoryImportService/import"
      interval="4000"
      name="DirectoryImportService"
      quarantine="quarantines/DirectoryImportService"
      root="roots/DirectoryImportService"/>
    <DicomImportService
      class="org.rsna.ctp.stdstages.DicomImportService"
      logConnections="no"
      name="DicomImportService"
      port="104"
      quarantine="quarantines/DicomImportService"
      root="roots/DicomImportService"/>
    <DirectoryStorageService
      acceptDuplicates="yes"
      class="org.rsna.ctp.stdstages.DirectoryStorageService"
      defaultString=""
      id="CollectionImportStorage"
      logDuplicates="no"
      name="DirectoryStorageService"
      quarantine="quarantines/DirectoryStorageService"
      root="roots/DirectoryStorageService"
      setStandardExtensions="no"
      structure="(0010,0020) - (0010,0010) / (0008,0020) / Series (0020,0011) "
      whitespaceReplacement="_"/>
  </Pipeline>
```

```

<Pipeline
  name="CollectionAnonymization"
  root="Collection/anonymization">
    <DirectoryImportService
      class="org.rsna.ctp.stdstages.DirectoryImportService"
      id="CollectionAnonymizerInput"
      import="roots/DirectoryImportService/import"
      interval="4000"
      name="DirectoryImportService"
      quarantine="quarantines/DirectoryImportService"
      root="roots/DirectoryImportService"/>
    <ObjectCache
      class="org.rsna.ctp.stdstages.ObjectCache"
      id="CollectionObjectCache"
      name="ObjectCache"
      root="roots/ObjectCache"/>
    <DicomAnonymizer
      class="org.rsna.ctp.stdstages.DicomAnonymizer"
      id="CollectionDicomAnonymizer"
      lookupTable="scripts/LookupTable.properties"
      name="DicomAnonymizer"
      quarantine="quarantines/DicomAnonymizer"
      root="roots/DicomAnonymizer"
      script="scripts/DicomAnonymizer.script"/>
    <ManifestLogger
      cacheID="CollectionObjectCache"
      class="edu.uams.tcia.ManifestLogger"
      id="CollectionManifestLogger"
      manifestLogID="CollectionManifestLog"
      name="ManifestLogger"
      root="roots/CollectionManifestLogger"/>
    <DirectoryStorageService
      acceptDuplicates="yes"
      cacheID="CollectionObjectCache"
      class="org.rsna.ctp.stdstages.DirectoryStorageService"
      defaultString=""
      id="CollectionAnonymizerStorage"
      logDuplicates="no"
      name="DirectoryStorageService"
      quarantine="quarantines/DirectoryStorageService"
      root="roots/DirectoryStorageService"
      setStandardExtensions="no"
      structure="(0010,0020)-(0010,0010)/(0008,0020)/Series(0020,0011)"
      whitespaceReplacement="_"/>
  </Pipeline>
  <Pipeline
    name="CollectionExport"
    root="Collection/export">
      <DirectoryImportService
        class="org.rsna.ctp.stdstages.DirectoryImportService"
        id="CollectionExportInput"
        import="roots/DirectoryImportService/import"
        interval="4000"
        name="DirectoryImportService"
        quarantine="quarantines/DirectoryImportService"
        root="roots/DirectoryImportService"/>
      <HttpExportService

```

```

        class="org.rsna.ctp.stdstages.HttpExportService"
        id="CollectionExportOutput"
        name="HttpExportService"
        quarantine="quarantines/HttpExportService"
        root="roots/HttpExportService"
        sendDigestHeader="no"
        url="http://192.168.0.225:7777"/>
    </Pipeline>
    <Pipeline
        name="DummyExportReceiver"
        root="Collection/dummy/dicom">
        <HttpImportService
            class="org.rsna.ctp.stdstages.HttpImportService"
            logConnections="no"
            name="HttpImportService"
            port="7777"
            quarantine="quarantines/HttpImportService"
            root="roots/HttpImportService"/>
        </Pipeline>
        <Pipeline
            name="DummyManifestReceiver"
            root="Collection/dummy/manifest">
            <HttpImportService
                class="org.rsna.ctp.stdstages.HttpImportService"
                logConnections="no"
                name="HttpImportService"
                port="7778"
                quarantine="quarantines/HttpImportService"
                root="roots/HttpImportService"/>
            </Pipeline>
        </Configuration>

```

The configuration has three main pipelines:

1. CollectionImport receives and stores DICOM files.
2. CollectionAnonymization anonymizes and stores files.
3. CollectionExport transmits files to the principal investigator.

The configuration has two additional pipelines just for testing:

1. DummyExportReceiver receives DICOM files from the CollectionExport pipeline.
2. DummyManifestReceiver receives manifest submissions.

The configuration has two Plugins:

- The TCIAPlugin provides the interface into the pipelines for the wizard.
- The ManifestLogPlugin accumulates manifest information for submission to the principal investigator.

3 User Accounts

When CTP starts, the TCIAPlugin changes the password of the **admin** user to **tcia** and grants it these roles:

- The **admin** role grants access to the administrative roles on the server.
- The **TCIA** role grants access to the TCIAServlet.
- The **qadmin** role grants access to the QuarantineServlet.
- The **shutdown** role allows the wizard to shut down CTP.

Additional users can be created through the User Manager on the CTP home page by logging in as a user with the **admin** role.

The wizard can automatically log in the user as tcia/tcia, or it can provide a UI that gets the username and password from the user and then make the call to the LoginServlet. The URL of the LoginServlet is:

/login/ajax?username=...&password=...

The login returns either a 200 or a 403 response code.

To provide a logout feature, the wizard can use this URL:

/login/ajax?logout

This call always returns 200.

4 The TCIAServlet API

This section describes the functions provided to the wizard by the TCIAServlet. All functions are accessed by making an HTTP connection to the CTP server. Unless otherwise indicated, all functions return an HTTP response with **Content-Type: text/xml; charset=UTF-8**.

4.1 Updating the DicomAnonymizer Lookup Table

To update the anonymizer lookup table from an Excel spreadsheet, the wizard does an HTTP POST of the file (with Content-Type **multipart/form-data**) to the URL:

/Collection

The spreadsheet must be an xlsx file with this structure:

	A	B	C	D
1		ptid	cname	ddate
2	Local Patient ID	Anonymized Patient ID	Collection Name	Diagnosis Date (M/D/YYYY)
3	657662	12345	A Collection Name	4/15/2014
4	5555	4444	Another Collection Name	2/1/2011
5	etc.			
6				

The first row contains the KeyTypes used in the lookup table for the data in each spreadsheet column. They must not be modified.

The second row provides a human-readable label for the column. The TCIA servlet needs the first row. It ignores the second row and starts processing rows at the third row.

4.2 Submitting Files to the Import Pipeline

4.2.1 Listing Files on the Server

To list the files in a directory on the server, the wizard does an HTTP GET to:

/Collection/listFiles?dir=path

The servlet returns an XML structure like this:

```
<dir
  name="...directory name..."
  parent="...absolute path to the parent directory..."
  <dir name="...child directory 1 name..." />
  <dir name="...child directory 2 name..." />
  ...
  <dir name="...child directory n name..." />
  <file name="...child file 1 name..." />
  <file name="...child file 2 name..." />
  ...
  <file name="...child file n name..." />
</dir>
```

4.2.2 Submitting Files to the DirectoryImportService

To submit a list of files to the DirectoryImportService in the import pipeline, the wizard does an HTTP GET to:

/Collection/submitFile?file="path sequence"

where **path sequence** is a string of paths separated by pipe characters, like this:

absolutePath-1|absolutePath-2|...|absolutePath-n

The servlet copies the files to the import pipeline. An absolutePath that points to a file imports that file. An absolutePath that points to a directory imports all the files in the directory and all of its child directories recursively.

The servlet returns an XML structure indicating whether the submissions all succeeded (<OK/>) or at least one failed (<NOTOK/>).

4.2.3 Submitting Files to the DicomImportService

The DicomImportService receives DICOM transfers on port 104. No configuration of the AE Titles is necessary; the SCP accepts all AE Titles.

4.3 Viewing the Status of the Import Pipeline

To obtain the number of files that have been received but not yet processed by the import pipeline, the wizard does an HTTP GET to the URL:

/Collection/getImportStatus

The AJAX call returns an XML structure like this:

```
<status queueSize="0" />
```

The queueSize attribute is the sum of the queue sizes of all the import services in the import pipeline. When the queueSize attribute is zero, all the received files have been passed down the pipeline and are ready for anonymization.

4.4 Viewing the List of Files Ready for Anonymization

To obtain the list of files that have been received but not yet anonymized, the wizard does an HTTP GET to the URL:

/Collection/listImport

The AJAX call returns an XML structure like this:

```
<?xml version="1.0" encoding="UTF-8"?>
<DicomFiles>
  <dir name="DirectoryStorageService" parent="D:\JavaPrograms\CTP-UAMS\CTP\Collection\import\roots">
    <dir name="1200824338-Bunny,Bugs">
      <dir name="20010312">
        <dir name="Series1">
          <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="1" StudyDate="20010312" name="9999.222149640329028447124359750781404748588"/>
          <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="1" StudyDate="20010312" name="9999.80349327100857023540751804102632703589"/>
        </dir>
      <dir name="Series2">
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.104115109181271829052136463305872175892"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.11385380932445884502388342839051934579"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.113983435978265066276300216885396006062"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.129003575003187817511384952750931955308"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.13364533319376454083949526458610875307"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.134095359512589080454306553825903370604"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.140285531251106151469597569011030160535"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.145415533035164672938705109621299252634"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.147559515087250199331639124618594546423"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.15693385102318300015440233716540985731"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.162760251557332918073608336500345336242"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.165810487962087107324577666065514793531"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.169831188606048016196033942139497761570"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.19145854833812192348065257738366545558"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.19451077797521877873061197566127973179"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.199259582581030369204776800763903847687"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.208689751371375329075102446940617839366"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.214257504674206793970526899765796655764"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.22094375288831707761401256395596647078"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.246112055202027104929794630877353288271"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.246505612395667702453661050403128845256"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.25925408076529146128857563062035589231"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.2670480403754384675701867718924305846358"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.273080104198635637751548367339060238477"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.2732522328372857128742264620878963063"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.28314757302252997677877678811457402258"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.28955223369947118261678031165462679034"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.289851792726025215340234985487577939151"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.292524086301371757006501187334377036996"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.30380247284068148351597079720694291413"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.36484771156478029866335762898287523006"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.5388163437211486952275333939442683389"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.72256014020957279464244896549265973992"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.73827376123490386449067540738220534566"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.7398455290608225677351279554331952611"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.78997342561141838539972752112481326238"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.88418205025288892784403384755226303426"/>
        <DicomObject Modality="CT" PatientID="1200824338" PatientName="Bunny,Bugs" Series="2" StudyDate="20010312" name="9999.95192804506073095562813808242229687492"/>
      </dir>
    </dir>
  </dir>
</DicomFiles>
```

The XML structure mimics the directory structure in which the files are stored. Each **dir** element in the XML structure represents a directory. The top-level directory is the root directory of the DirectoryStorageService pipeline stage that contains the files received by the import pipeline. Under the top-level directory, there is one directory for each patient. Under a patient's directory is one directory for each date on which the patient had a study. The study date directory contains all the images for studies done on that date. Individual image files are represented in the XML structure by **DicomObject** elements.

The relative path to a directory can be obtained by walking the tree from the point in question back to the top. Thus, the relative path to the **20010312** directory is:

DirectoryStorageService/1200824338-Bunny,Bugs/20010312

The base of the relative path is the root of the pipeline stage.

4.5 Submitting Files for Anonymization

To anonymize studies that were listed in 4.2, above, the wizard passes some level in the hierarchy (a patient, a single study, or the top-level directory) to the servlet in an HTTP GET to the URL:

/Collection/anonymize?file=filepath

where **filepath** is the relative path. If the filepath is a directory, the servlet processes all the files in that directory and all its child directories. For example, to process everything that has been received, the URL would be:

/Collection/anonymize?file=DirectoryStorageService

The servlet moves the files from the import pipeline to the anonymizer pipeline to start the anonymization. The servlet returns an XML structure with one element <OK/> or <NOTOK/> to indicate whether all the moves were successful.

4.6 Viewing the List of Files Ready for Export

To obtain the list of files that have been anonymized but not yet exported, the wizard does an HTTP GET to the URL:

/Collection/listAnonymized

The AJAX call returns an XML structure like this:

```
<?xml version="1.0"?>
<DicomFiles>
  <dir name="/DirectoryStorageService" parent="D:\JavaPrograms\CTP-UAMS\CTP\Collection\anonymization\roots">
    <dir name="1200824338-Bunny_Bugs">
      <dir name="20010312">
        <dir name="Series1">
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="1" StudyDate="19920528" name="1.2.840.113654.2.70.1.196154021349712661720723242598417253310"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="1" StudyDate="19920528" name="1.2.840.113654.2.70.1.24987608175328106310133805613683356439"/>
        </dir>
        <dir name="Series2">
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.100668460971350414894490102932760453237"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.10206093558354332079943683094877159893"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.10577958986487253224308501378982492719"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.11549220012628587817540242011917783060"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.11717510016388387727276475378212702737"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.123921539284844847827383058646932152096"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.128675395687132861321061084828296976924"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.130421976200204254919547971728858992611"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.13944181554570888557228572515159409240"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.14258958814260588904253257014007813483"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.14594928990519624317297006636062544279"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.159422924702382504851306269862036196064"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.161158294508295584409728267395199962999"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.161699679064174416207625600477645061173"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.16823997610907354322962188414701280751"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.187189977615953577403173506995370930482"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.198466521145949332279162363066595489659"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.207220241178249664138652977244491393663"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.246463463499649159607524252287593903016"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.24658802765830245884215749598059574558"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.259485621065070668421580885329452131524"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.26257742661187349411641009557696652205"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.267271874089169982177365211904539741389"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.267680436852985977609972671988518120083"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.269659613708517497492872534597899020830"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.2773805393474173938308147833420860646925"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.305943958418409172941668374556534215576"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.31369817994981755287724376832201741454"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.32627568491443177480275684640795946119"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.3267035388799401327085691634463164481"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.60069941867498313588041854603472491336"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.6167432526989003066455558088229901963"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.667619360275242009041302999923054336"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.68829204219636694715410368918017413798"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.69454031442707587359104268938477239058"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.696871331479752617938799266983290874"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.774554247789403558338619291800497694"/>
          <DicomObject Modality="CT" PatientID="3654257458" PatientName="204173" Series="2" StudyDate="19920528" name="1.2.840.113654.2.70.1.99212405854561686185313042515920356"/>
        </dir>
      </dir>
    </dir>
  </DicomFiles>
```

The XML structure exactly mimics the one that listed the imported files, but it references a different DirectoryStorageService pipeline stage, so although the directory names are the same, they are relative to a different root directory.

[This is probably also a good place to provide file viewing capabilities to reassure the user that they are not exporting PHI.]

4.7 Exporting Files

To export studies that were listed in 4.4, above, the wizard passes some level in the hierarchy (a patient, a single study, or the top-level directory) to the servlet in an HTTP GET to the URL:

/Collection/export?file=filepath

where **filepath** is again the relative path.

The servlet moves the files from the anonymizer pipeline to the export pipeline to start the export. The servlet returns an XML structure with one element <OK/> or <NOTOK/> to indicate whether all the moves were successful.

4.8 The Manifest

As files are anonymized, a manifest listing key identifiers from each file is created. The manifest can be obtained as either XML or CSV, with or without PHI. The manifest contains one entry for each series. The TCIAServlet provides five functions for accessing the manifest:

4.8.1 Clearing the Manifest

The wizard can clear the manifest by an HTTP GET to the URL:

/Collection/clearManifest

The servlet returns an XML structure with one element <OK/> or <NOTOK/> to indicate whether all the operation was successful.

4.8.2 Getting the Manifest Status

The wizard can get an XML structure containing the status of the processing of files in the Anonymizer pipeline by an HTTP GET to the URL:

/Collection/getManifestStatus

The servlet returns an XML structure like this:

```
<Status
  startingQuarantineCount="0"
  currentQuarantineCount="0"
  queuedInstanceCount="40"
  currentManifestInstanceCount="40"/>
```

4.8.3 Listing the Manifest with PHI as XML

The wizard can obtain an XML structure containing the manifest of series, including PHI, by an HTTP GET to:

/Collection/listLocalManifest/xml

The resulting structure looks like this:

```

▼<Manifest>
  ▼<Series>
    <Collection value=""/>
    <SiteName value=""/>
    <PatientID phi="1200824338" value="3654257458"/>
    <SOPClassUID value="1.2.840.10008.5.1.4.1.1.2"/>
    <Modality value="CT"/>
    <StudyDate phi="20010312" value="19920528"/>
    <SeriesDate phi="" value=""/>
    <StudyDescription value=""/>
    <SeriesDescription value="Sinus Coronal"/>
    <StudyInstanceUID phi="9999.294490511222623278153251388605292076811" value="1.2.840.113654.2.70.1.178587914247804709210381252427058545955"/>
    <SeriesInstanceUID phi="9999.13474336704961777217271177118548448960" value="1.2.840.113654.2.70.1.113479017276406998463382173796105434129"/>
    <NumFiles value="38"/>
  </Series>
  ▼<Series>
    <Collection value=""/>
    <SiteName value=""/>
    <PatientID phi="1200824338" value="3654257458"/>
    <SOPClassUID value="1.2.840.10008.5.1.4.1.1.2"/>
    <Modality value="CT"/>
    <StudyDate phi="20010312" value="19920528"/>
    <SeriesDate phi="" value=""/>
    <StudyDescription value=""/>
    <SeriesDescription value="Topogram 0.6 T20s"/>
    <StudyInstanceUID phi="9999.294490511222623278153251388605292076811" value="1.2.840.113654.2.70.1.178587914247804709210381252427058545955"/>
    <SeriesInstanceUID phi="9999.264302660968596814195213146253418668559" value="1.2.840.113654.2.70.1.116933866301694768395323515093764104644"/>
    <NumFiles value="2"/>
  </Series>
</Manifest>

```

4.8.4 Listing the Manifest without PHI as XML

The wizard can obtain an XML structure containing the manifest of series, including PHI, by an HTTP GET to:

/Collection/listExportManifest/xml

The resulting structure looks like this note that the phi attributes of the elements are not present):

```

▼<Manifest>
  ▼<Series>
    <Collection value=""/>
    <SiteName value=""/>
    <PatientID value="3654257458"/>
    <SOPClassUID value="1.2.840.10008.5.1.4.1.1.2"/>
    <Modality value="CT"/>
    <StudyDate value="19920528"/>
    <SeriesDate value=""/>
    <StudyDescription value=""/>
    <SeriesDescription value="Sinus Coronal"/>
    <StudyInstanceUID value="1.2.840.113654.2.70.1.178587914247804709210381252427058545955"/>
    <SeriesInstanceUID value="1.2.840.113654.2.70.1.113479017276406998463382173796105434129"/>
    <NumFiles value="38"/>
  </Series>
  ▼<Series>
    <Collection value=""/>
    <SiteName value=""/>
    <PatientID value="3654257458"/>
    <SOPClassUID value="1.2.840.10008.5.1.4.1.1.2"/>
    <Modality value="CT"/>
    <StudyDate value="19920528"/>
    <SeriesDate value=""/>
    <StudyDescription value=""/>
    <SeriesDescription value="Topogram 0.6 T20s"/>
    <StudyInstanceUID value="1.2.840.113654.2.70.1.178587914247804709210381252427058545955"/>
    <SeriesInstanceUID value="1.2.840.113654.2.70.1.116933866301694768395323515093764104644"/>
    <NumFiles value="2"/>
  </Series>
</Manifest>

```

4.8.5 Listing the Manifest with PHI as CSV

The wizard can obtain a text string containing the manifest as a spreadsheet in CSV format by an HTTP GET to:

/Collection/listLocalManifest/csv

This AJAX call returns a CSV text response with **Content-Type: text/csv; charset=UTF-8**. The resulting spreadsheet has two rows for each series, one for the anonymized values and one for the PHI values. Series are separated in the spreadsheet by blank lines.

4.8.6 Listing the Manifest without PHI as CSV

The wizard can obtain a text string containing the manifest as a spreadsheet in CSV format by an HTTP GET to:

/Collection/listExportManifest/csv

This AJAX call returns a CSV text response with **Content-Type: text/csv; charset=UTF-8**. The resulting spreadsheet contains one row for each series, with only the anonymized values.

4.8.7 Listing the Manifest with PHI as XLSX

The wizard can obtain a text string containing the manifest as a spreadsheet in XLSX format by an HTTP GET to:

/Collection/listLocalManifest/xlsx

This AJAX call returns a CSV text response with **Content-Type: application/vnd.openxmlformats-officedocument.spreadsheetml.sheet**. The resulting spreadsheet has two rows for each series, one for the anonymized values and one for the PHI values. Series are separated in the spreadsheet by blank lines.

4.8.8 Listing the Manifest without PHI as XLSX

The wizard can obtain a text string containing the manifest as a spreadsheet in XLSX format by an HTTP GET to:

/Collection/listExportManifest/xlsx

This AJAX call returns a CSV text response with **Content-Type: application/vnd.openxmlformats-officedocument.spreadsheetml.sheet**. The resulting spreadsheet contains one row for each series, with only the anonymized values.

4.8.9 Exporting the Manifest

The wizard can send the CSV manifest without PHI to the principal investigator by an HTTP GET to:

/Collection/exportManifest

The servlet returns an XML structure with one element <OK/> or <NOTOK/> to indicate whether the transmission was successful.

4.9 Image Functions

4.9.1 Listing the Elements of a DICOM File

The wizard can obtain an HTML page containing a table showing the values of the elements in a DICOM file by an HTTP GET to the URL:

/Collection/listElements?file=filepath

The servlet returns a JPEG image in the response with **Content-Type: image/jpeg**. If the servlet cannot find the file, the servlet returns a 404 response code. If the servlet cannot parse the file, the servlet returns a 500 response code and enters a stack trace in the log.

4.9.2 Getting a JPEG Image from a DICOM File

The wizard can obtain a browser-viewable image of a DICOM by an HTTP GET to the URL:

/Collection/getImage?file=filepath

The servlet returns a JPEG image in the response with **Content-Type: image/jpeg**. If the servlet cannot find the file, the servlet returns a 404 response code. If the servlet cannot parse the file, the servlet returns a 500 response code and enters a stack trace in the log. If a 500 is received and the log shows a **No class def found** error for the StreamSegmentMapper, it means that either the ImageIO Tools have not been installed, or the image has a format is not supported by the ImageIO Tools.

4.10 Special Functions

4.10.1 Getting the List of File System Roots on the Server

The wizard can obtain an XML structure listing the file system roots on the server by an HTTP GET to the URL:

/Collection/getFileSystemRoots

This servlet returns an XML structure like this:

```
<roots>
  <root name="C:\" desc="Local Disk"/>
  <root name="D:\" desc="Local Disk"/>
  <root name="E:\" desc="Local Disk"/>
  <root name="P:\" desc="CD Drive"/>
  <root name="Q:\" desc="CD Drive"/>
</roots>
```

4.10.2 Getting the Available Space on the Server

The wizard can obtain an XML structure indicating the available space on the partition on which the CTP instance is located by an HTTP GET to the URL:

/Collection/getAvailableSpace

This servlet returns an XML structure like this:

```
<space partition="D:\" available="434932" units="MB"/>
```

The wizard can specify the file system root on which to obtain the available space by including the **root** query parameter:

/Collection/getAvailableSpace?root=...

Where the value of the **root** parameter is a name obtained from the `getFileSystemRoots` function.

4.10.3 Getting the Export Queue Size

The wizard can obtain an XML structure indicating the current size of the export queue by an HTTP GET to the URL:

/Collection/getExportQueueSize

This servlet returns an XML structure like this:

```
<queue stage="HttpExportService" size="10"/>
```

4.10.4 Getting the URL of the Quarantine Manager Servlet

The wizard can obtain an XML structure containing the URL of the CTP Quarantine Manager for the DicomAnonymizer quarantine by an HTTP GET to the URL:

/Collection/getQuarantineURL

This servlet returns an XML structure like this:

```
<quarantine stage="DicomAnonymizer" url="/quarantines?p=1&s=2"/>
```

4.10.5 Getting a Summary of the DicomAnonymizer Quarantine

The wizard can obtain an XML structure containing a summary of the objects in the DicomAnonymizer quarantine by an HTTP GET to the URL:

/Collection/getQuarantineSummary

This call returns an XML structure like this:

... TBD ...

4.10.6 Shutting Down the Server

The wizard can shut down the CTP server by an HTTP GET to:

/shutdown

The system returns either a 200 or 403 response code. It fails if the user does not have the **shutdown** role. This URL accesses the CTP ShutdownServlet directly. It returns a web page like this:

```
Shutdown request received from admin at 192.168.0.225.  
Goodbye.
```

4.11 Testing Functions

During testing, it may be convenient to clear all the import/export directories, the DicomAnonymizer quarantine, and the manifest. This can be done by an HTTP GET to the URL:

/Collection/reset

The servlet does all it can and returns an XML structure with one element <OK/>. It never returns <NOTOK/>.